


ALL SUTURE SOFT ANCHOR

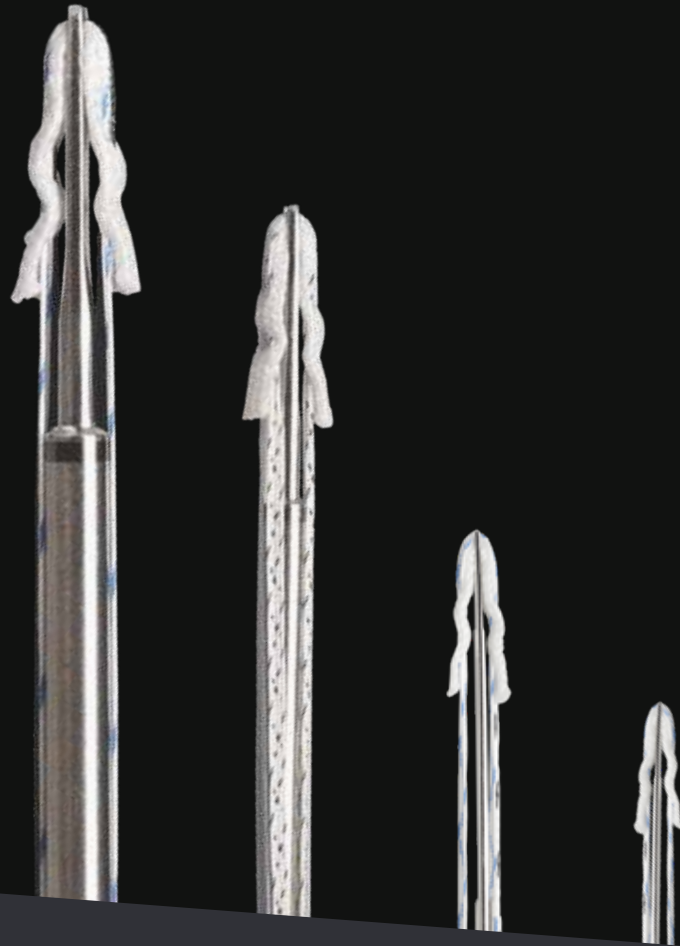

ALL SUTURE ANCHOR




ALL SUTURE ANCHOR

It is specially designed for the fixation of soft or prosthetic tissues to the bone

SOFT ANCHOR are impacted anchor of small diameter exclusively composed of textile fibers.

The polyester body of the anchor is associated with one or two UHMWPE suture threads which are mounted on a disposable screwdriver allowing insertion into the bone.

The traction on the suture threads makes it possible to modify the shape of the body of the anchor which thus locks itself behind the cortical bone and ensures an effective anchorage.

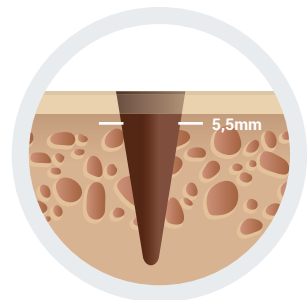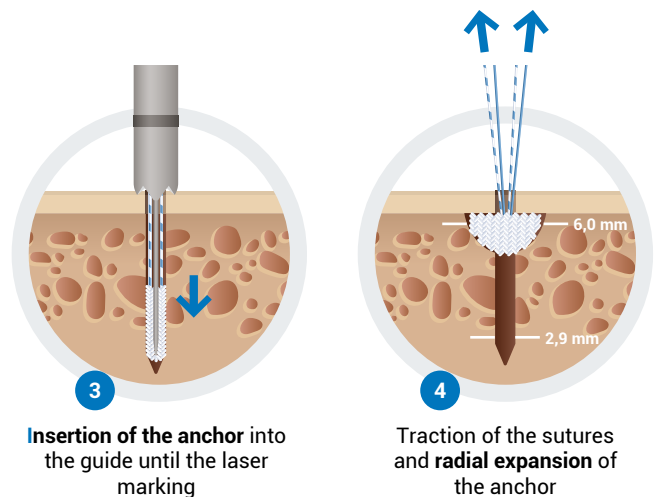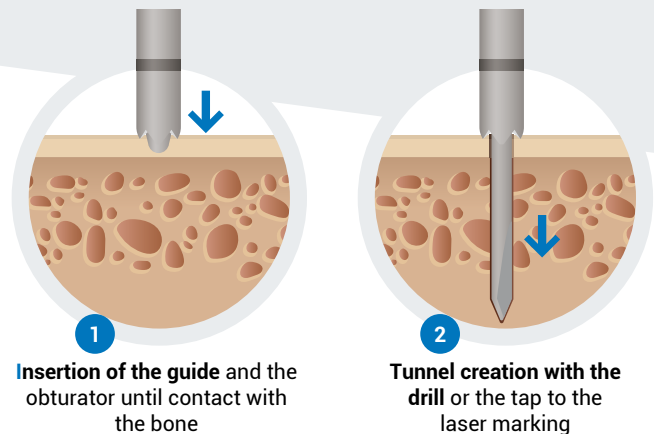

Bone economy compared to a screwed anchor

ADVANTAGES

- Bone economy
- Easy to use
- No rigid material in case of revision or in the joint in case of tearing
- Radio-transparent
- Can be used under arthroscopy or mini open (with needles)
- Biocompatible
- High pull-out resistance
- Non-absorbable

Beveled tip for better impaction in bone tissue



ALL SUTURE ANCHOR



SHOULDER
ROTATOR CUFF,
BICEPS TENODESIS AND
REINSERTION

Ø2,9
mm



SHOULDER
SLAP, BANKART
INSTABILITY



ANKLE
TENDON AND LIGAMENT
REINSERTION

Ø1,8
mm



HAND
TENDON AND LIGAMENT
REINSERTION

Ø1,0
mm



ANKLE
TENDON AND LIGAMENT
REINSERTION

Ø1,5
mm



SUTURE ADVANTAGES

ARTHROVIMS® SOFT ANCHOR are associated with 1 or 2 sutures, sutures made of very high molecular weight polyethylene (UHMWPE)

HIGH TENSILE STRENGTH
EASY REALIZATION OF KNOTS
STRENGTH AND SAFETY OF SUTURES



ALL SUTURE ANCHOR

STORAGE

SOFT ANCHOR must be stored in its original unopened packaging, in a clean and dry place, at ambient temperature.

PRODUCTS & INSTRUMENTATION

SOFT ANCHOR

is packaged in a double pouch and sterilized by ethylene oxide. Instrumentation for arthroscopic procedures is available. Single use. **Do not re-sterilize.**

SOFT ANCHOR must be used with the corresponding instrumentation supplied by NEOSYS. For any further information, please refer to the IFU.

NeoSys

Medical Solutions

WWW.NEOSYSMED.COM

INSTRUMENTATION

SOFT ANCHOR must be used with the specific instrumentation provided by NEOSYS.

CE 0653
CLASS IIb / ANCHOR

CE
CLASSE I / INSTRUMENTATION

